

Create your study template on the WOMBAT web application



Users

[Sign Up](#) | [Login](#) | [Forgot your password?](#)

Please enter your email and password to log in


Email

Password

Login

Step 4a: Create a new study

You are logged in as [Edit My Details](#) | [Change My Password](#) | [Logout](#)



Administration Manage Users	Study Configuration List Create Import	Study Data Download
---------------------------------------	---	-------------------------------

Create New Study

Name

Session duration (minutes)

End session alert?

Step 4b: Add your dimensions

Configure

Session duration

Please click

Add Dimension

Add Dimension

Name

Multiple Selections

Mandatory

Save Cancel

Dimensions

Limitations

You can have as many dimensions, categories and subcategories as you like. However once you reach a certain number, the user will have to scroll on the tablet to view all the options.

Below are the approximate limits:

- With 4 dimensions, you can fit approximately 30 categories on screen at one time (spread across the 4 dimensions)
- Category names are truncated at 12 characters on the tablet
- You can fit approximately 9 subcategories on the screen at one time

Step 4c: Add your categories

What Edit Delete Please note the categories and subcategories in this dimension will be used as labels in the task list on

Single/Mu
Mandatory
Default Ca
Add Catego

Limitation

You can hav
number, the

Below are th

Edit Categories

Direct Care	<u>Delete</u>
Indirect Care	<u>Delete</u>
Medication	<u>Delete</u>
Documentation	<u>Delete</u>
Professional Communication	<u>Delete</u>
Administration	<u>Delete</u>
In transit	<u>Delete</u>
Supervision/Education	<u>Delete</u>
Social	<u>Delete</u>

Add more

Save Cancel

- dimensions)
- Category names are truncated at 12 characters on the tablet
- You can fit approximately 8 subcategories on the screen at one time

Step 4c: Add any subcategories

What [Edit](#) [Delete](#) *Please note the categories and subcategories in this dimension will be used as labels in the task list on the tablet.*

Single/Multiple Selection(s): Single

Mandatory/Optional: Mandatory

Default Category: None

[Edit Categories](#)

Direct Care [Add Subcategories](#) [Make default](#)

Indirect Care [Add Subcategories](#) [Make default](#)

Medication [Add Subcategories](#) [Make default](#)

Documentation [Add Subcategories](#) [Make default](#)


Professional Communication [Add Subcategories](#) [Make default](#)

Administration [Add Subcategories](#) [Make default](#)

In transit [Add Subcategories](#) [Make default](#)

Step 4d: Finish designing study template

You are logged in as [Edit My Details](#) | [Change My Password](#) | [Logout](#)



Administration Manage Users	Study Configuration List Create Import	Study Data Download
---------------------------------------	--	-------------------------------

The subcategories were successfully updated.

Configuration for the Doctors Work Measurement V1 study

Session duration: 2 hours

[Add Dimension](#)

[Finished](#)

Sort Dimensions


What [Edit](#) [Delete](#) Please note the categories and subcategories in this dimension will be used as labels in the task list on the tablet.



Done

Step 4e: Mark template as 'Active'

You are logged in as [Edit My Details](#) | [Change My Password](#) | [Logout](#)



Administration Manage Users	Study Configuration List Create Import	Study Data Download
---------------------------------------	--	-------------------------------

Manage Study Doctors Work Measurement V1

Description

Mark the status of your template, then click 'Save'



- Exemplar
- Active
- Archived

Done




Downloading a study template

Click on the
WOMBAT

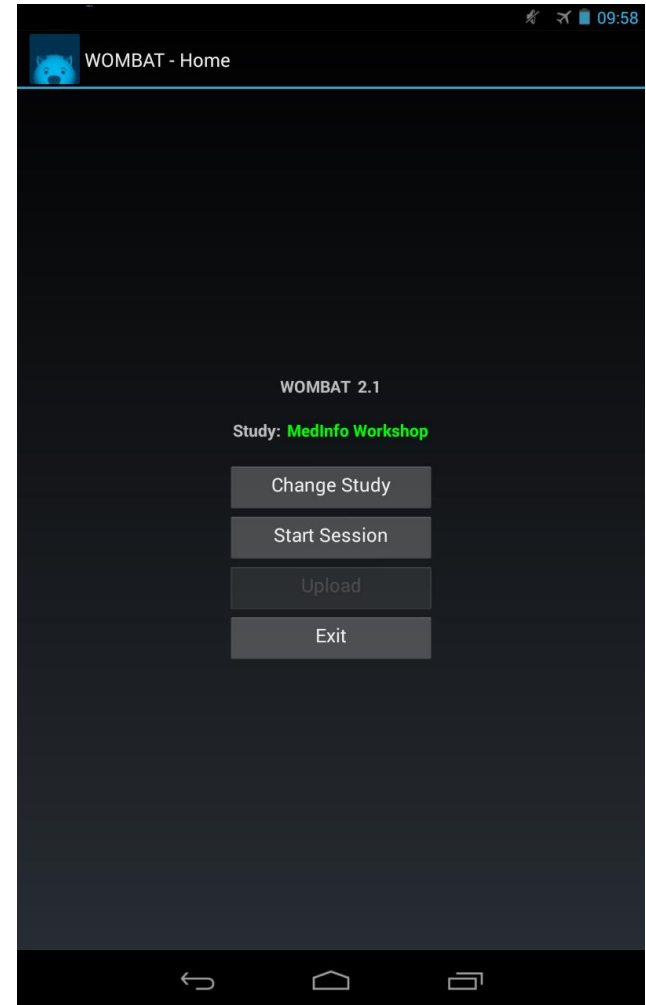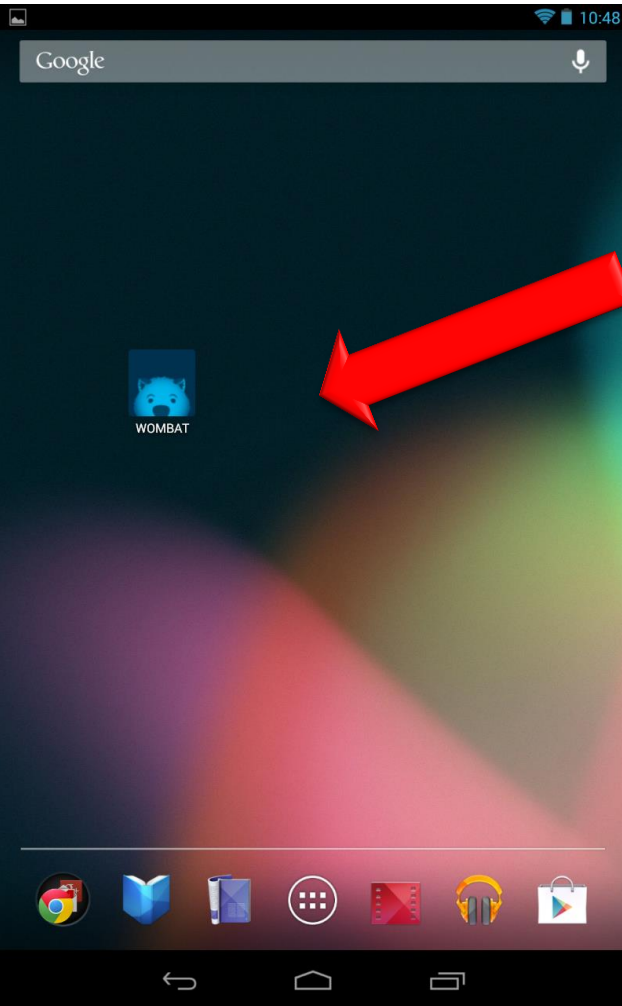

Starting a session

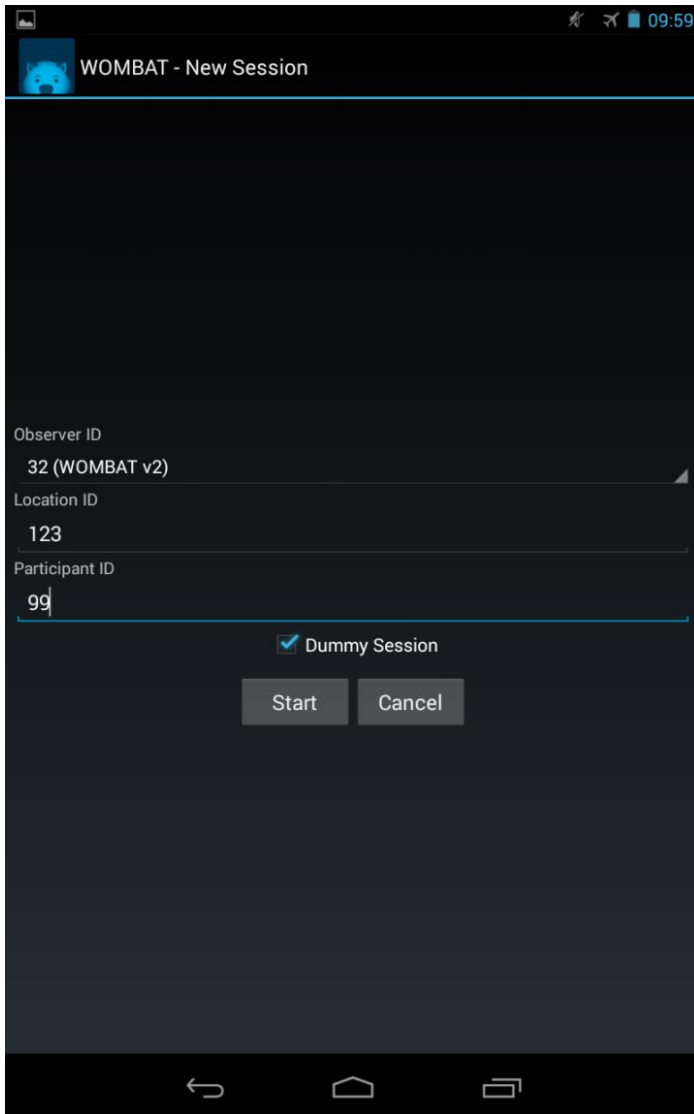

The screenshot shows the 'WOMBAT - New Session' screen on a mobile device. The status bar at the top indicates the time is 09:59. The app title 'WOMBAT - New Session' is displayed in the top left. Below the title, there are three input fields: 'Observer ID' with the value '32 (WOMBAT v2)', 'Location ID' with the value '123', and 'Participant ID' with the value '99'. A checkbox labeled 'Dummy Session' is checked. At the bottom of the form, there are two buttons: 'Start' and 'Cancel'. The Android navigation bar is visible at the very bottom.

Timer starts
when you press
“start”

